



# Monofocal IOLs

## Aspheric HP IOLs

\* SeeLens HP

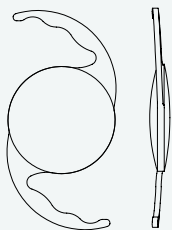

\* BunnyLens HP

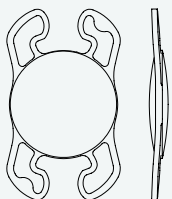

\* Available in preloaded configuration

	SeeLens HP *	BunnyLens HP *
<b>Platform</b>	C-Loop	4 - Loop
<b>Overall Diameter</b>	13 mm	11 mm (>10D) 11.5mm (≤10D)
<b>Placement</b>	Capsular Bag	
<b>Optic diameter</b>	6 mm	
<b>Power range</b>	+5.0 to +10.0 (1D increments) +10.5 to +30.0 (0.5D increments) +31.0 to +35.0 (1D increments)	
<b>Optic design</b>	Aspheric	
<b>360° Continuous Square Edge</b>	Yes	
<b>Haptic angulation</b>	5°	
<b>Material</b>	Hydrophobic acrylic with bonded UV absorber and violet light filter	
<b>Refractive Index</b>	1.48 (@ 35°C)	
<b>A-constant (SRK/T) for Optical or Immersion US biometry</b>	119	118.9
<b>A-constant (SRK/T) for Contact US biometry</b>	118.5	118.4
<b>Sterilization</b>	ETO	


# Providing Excellent Vision Quality For Your Patient

The Aspheric HP IOLs were designed by the most advanced optical tools to provide:

- reduced spherical aberration
- functional vision in night conditions

The Aspheric Hydrophobic lens from Hanita Lenses, provides the patient with an excellent vision quality at day and night conditions, by using state-of-the-art aberration free Aspherical optical design.

## Stability and Centration

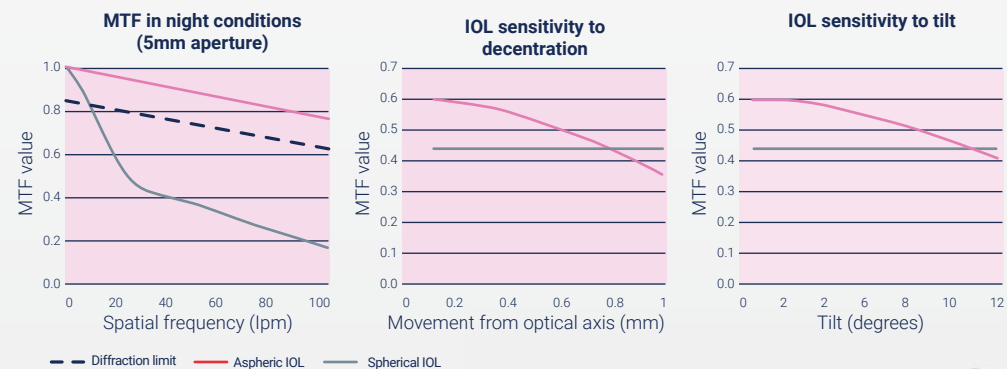

## PRELOADED HYDROPHOBIC IOL SYSTEM

- Simplifying cataract surgery
- Ready for implant
- Smooth and EASY injection through a 2.2 mm incision

